


OPT-DPP-COIN-S

Off-line processing by software
 Coincidence processing is performed on the output list data.

Repeated analysis is possible using list data.

- Example of use**
- #1 Specify the read data size.
 - #2 Specify coincidence time.
 *The unit is digit, 1digit = 3.90625ps.
 - #3 Specify the target channel for coincidence.

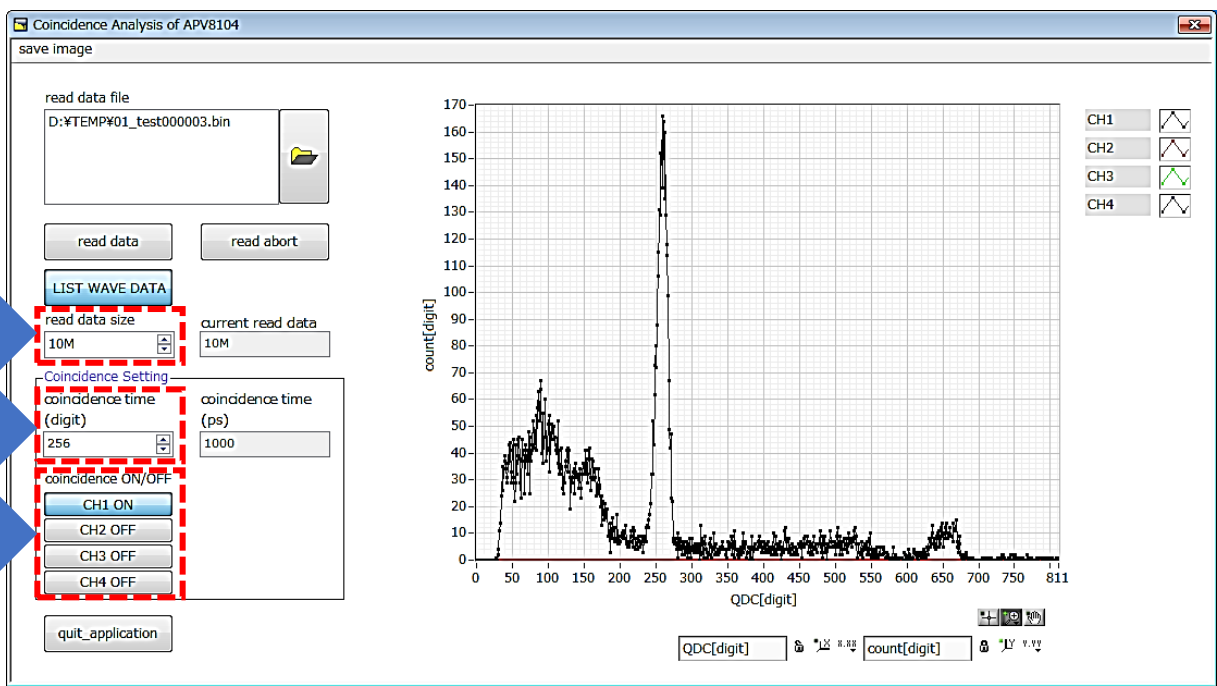

Mountable DPP

APV8104 APV8508	APV8108 APV8516	APV82516	APV8M24 APV8M44
OK	OK	OK	OK

Customizable upon customer request

*Images is for illustration purpose.
 *Please note that contents may change without prior notice.

